# **geminipro**

ITALY

## **20-IN-1 WIRELESS SMART BODY FAT SCALE**

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**MODEL NO 型號: GFS90S**

Before using your product for first time,  
please read this instruction manual.  
在使用本產品前，請仔細閱讀本說明書並保存好。

# IMPORTANT SAFEGUARDS

- **READ ALL INSTRUCTIONS BEFORE USE**

- **Pregnant woman and users with a pacemaker or electric equipment inside the body shall not use the scale.**
- The glass surface of the scale should be kept dry to prevent slipping.
- The scale can only bear normal static pressure when it is weighing. Under abnormal use or impact, there may be risk of damage to the scale.
- This scale is not waterproof; never immerse the scale in water or other liquid.
- Clean the scale with a damp cloth but do not allow water to go inside the scale and do not use any chemical cleaning agents.
- Remove the batteries if the scale is not being used for a long time. Never place any heavy objects on the equipment.
- The scale is for family use to measure weight, body fat, water mass, muscle mass, bone mass, calorie, BMI. It is not used for commercial or other weighing purpose.
- The body fat parameter is only for your reference. If you target to loss or gain weight, consult your doctor for medical advice.
- In order to get the accurate parameter, please make sure the surrounding area of scale has no electromagnetic objects to avoid interference.
- Please do not repair this scale by yourself. If the scale does not operate normally, check the batteries condition and if the connection direction is right. If the problem cannot be solved, please contact our after-sales service center for assistance.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- The minimum body scale starts at 5kg, which is not suitable for weighing objects smaller than 5kg.

# PARTS IDENTIFICATION

Before using the scale for the first time, please be familiar with all of the parts. Read all Instructions and Safeguards carefully.

LCD Display

BACK PANEL

Unit button


## 4 electrode on the glass platform

### The correct Fat/water mass measurement method:

Please stand on the scale barefooted, keep your soles of the feet clean and dry and has a good contact with the precision sensors.


Please refer to the following operation steps for the result.
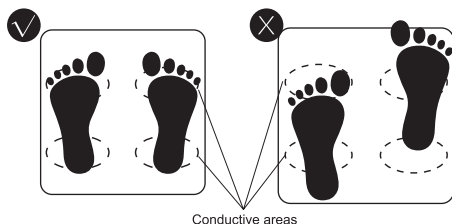
### Battery slot and cover

Lift the battery cover at the bottom of the scale and insert 3 x AAA batteries.

## LCD Display content

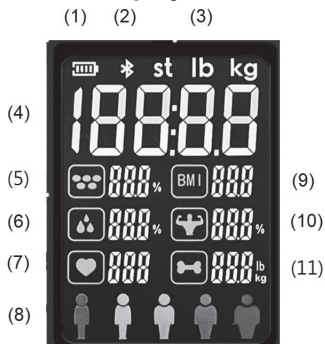- (1) Power level
- (2) Bluetooth connection
- (3) Unit – st / lb / kg
- (4) Weight
- (5) Fat mass
- (6) Body water
- (7) Heart rate
- (8) Obesity level
- (9) BMI
- (10) Muscle Rate
- (11) Bone Mass

## AUTOMATIC SHUT-OFF

To conserve battery life this scale has an automatic shut-off feature. If the scale is turned on and not operating, it will display “0.0” and turn off after 20 seconds. When the scale is turned on, a measurement is taken, and weight is locked, other measurement data will be displayed and the scale will shut off after 20 seconds when all datas are display completely.

## USING THE APPLIANCE

- Before using the scale, make sure you have entered all your personal data through the Leaone app.
- Put the scale on a hard and flat floor, avoid to use it on carpet or other soft surface, please gently step on the scale to wake it up and wait at least 10 seconds before weighing on the scale. You must also wake up the scale if it has been moved. Otherwise it will influence the accuracy.
- Use the scale while barefoot. The scale cannot measure body fat percentage if you're wearing shoes or socks.
- Stand on the scale directly when using the weight function only.
- To select the measuring unit at the beginning, press the “Unit” button which is located at the scale's base, the corresponding unit can be selected on the LCD screen, i.e. Kilograms (kg), Pounds (lb) or Stone: Pound (st:lb). Choose the desired “Unit” to start the measurement.
- 1 Stone(st)=14 Pounds(lb), 1 Kilogram(kg)=2.2 Pounds(lb). If a person's weight is displayed as 7:12st when using the st unit, the weight it represents is 7st+12lb
- Each time it is enabled, the first data is the Wake-up function, the second and third data shall prevail.

## MOBILE DEVICE SUPPORT

iOS 13.0 and above

Android 6.0 and above with Bluetooth 4.0

## APP INSTALLATION


1. Scan the QR Code or download the Leaone app from the Apple Store or Google Play.
2. Download and install the Leaone app on your device.
3. The App icon will appear on your phone or tablet after the installation is completed.

# INSTRUCTIONS FOR USE

## Step 1

Open the Leaone app on your device  
\*Please allow to search for Bluetooth on your devices


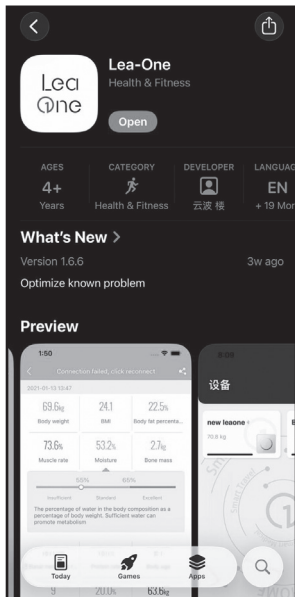
## Step 2

Register your account by email or login with social media account.


1. Click "Register" to registration page
2. Enter email address
3. Input password
4. Press "Register"

*P.S. This explanation is under the iOS version, for the android version, please follow the instruction on the app.*
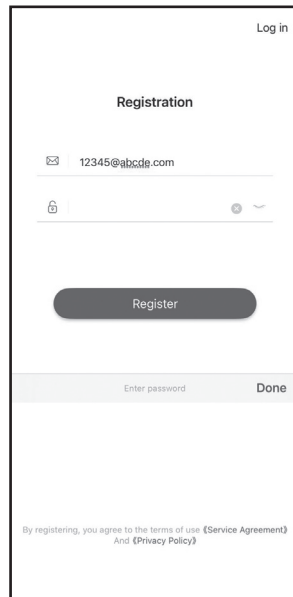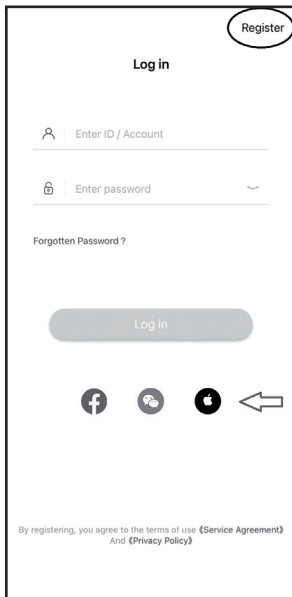
## Step 3

Fill in your personal details and press "OK"

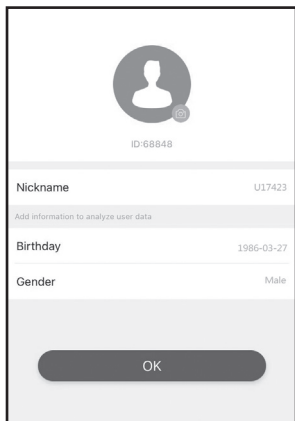**Add device – Two ways for connecting**

1) Bluetooth searching the devices nearby (Please remember that the Bluetooth and the scale must keep to wake up for connection on the processing)

**Step 1**

Press “Add device” for searching

**Step 2**

The device will shown on the screen once searched and press the icon for connect.

**Step 3**

Connecting


**Step 4**

Once the App connected with the device, it will show


**Step 5**

Press “Start using” once connected

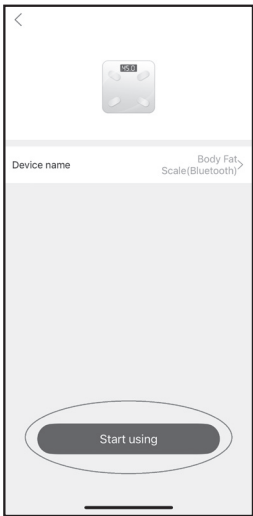**Step 6**


Choose unit and press “OK” to start for measurement.

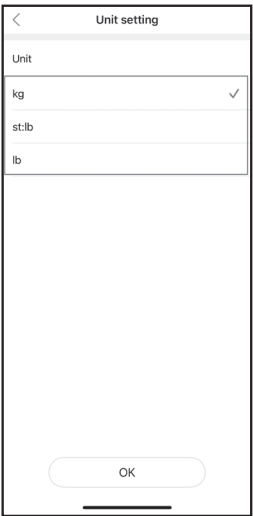2) Add manually (Please remember that the Bluetooth and the scale must keep to wake up for connection on the processing)


**Step 1**

Press “Add device” for searching

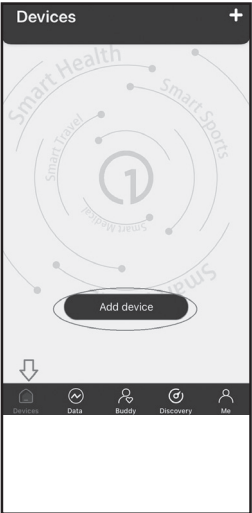**Step 2**

Press “Add manually” and then choose the device “Body fat scale (Bluetooth)”

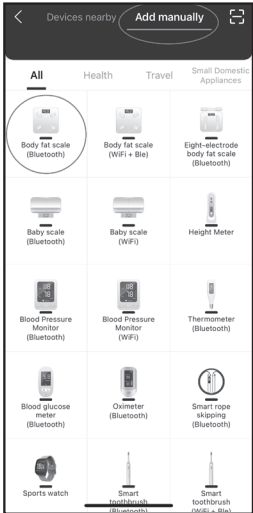
**Step 3**

Searching and you will see the Body Fat Scale (Bluetooth) on the screen and press it for connection.


**Step 4**

Connecting and you will see matched successfully once device is connected

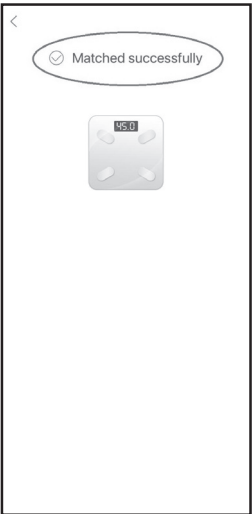
**Step 5**

Press “Start using” once connected

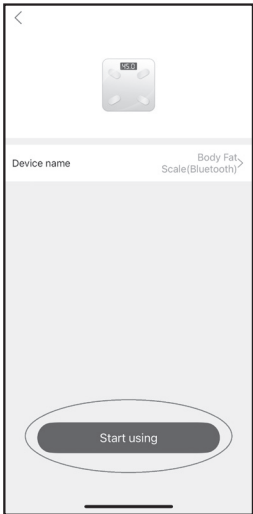**Step 6**


Choose unit and press “OK” to start for measurement.
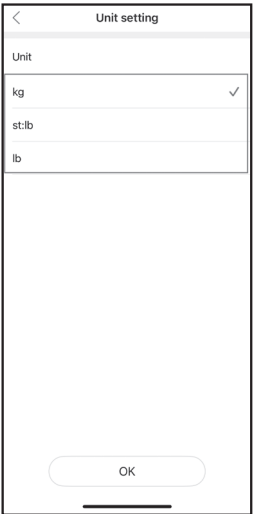
# TEST 20 BODY DATA IN 1 MEASUREMENT

Test 20 body data: Weight/ BMI/ BFR/ Muscle Rate/ Body water/ Bone mass/ BMR/Protein rate/ Metabolic age/ Visceral fat index/ Subcutaneous fat/ Standard weight/Weight control/ Fat mass/ Weight without fat/ Muscle mass/ Protein mass/ Obesity level/Body type/ Heart rate

The connection status will be shown “Connected” once the device connection successful. Then you can start to measurement.  
(PS: only 1 device can be connected at the same time.)




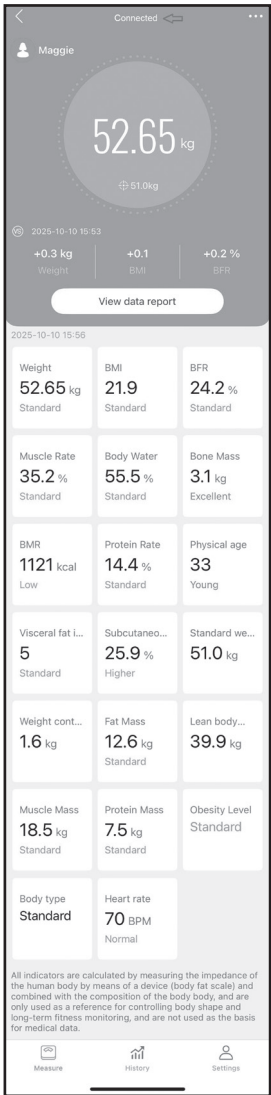
# EVALUATION OF RESULTS

After client the measuring data, the status will show “Measurement completed” and it will show the evaluation of results in the App or export your body composition analysis report to PDF at “View data report”

Press “history” icon to view your history data

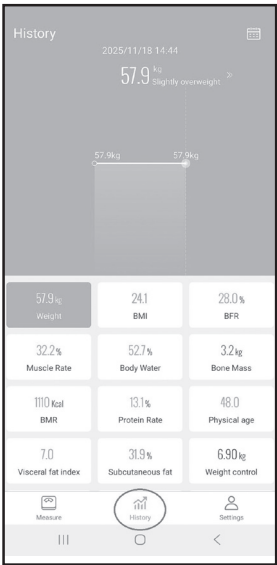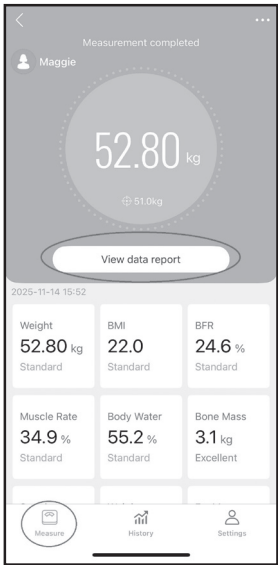
Press “Calendar” icon to view your data on specific date
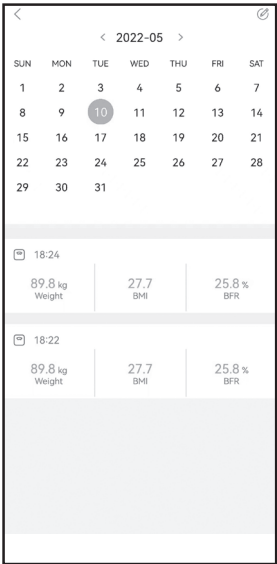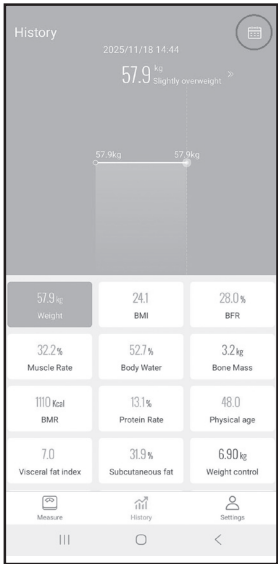
## OTHER FUNCTIONS

3 functions on the “settings” icon

- 1) About device – for setting unit or remove the device
- 2) Themes - set the theme color
- 3) Measurement reminder – set the time for reminder

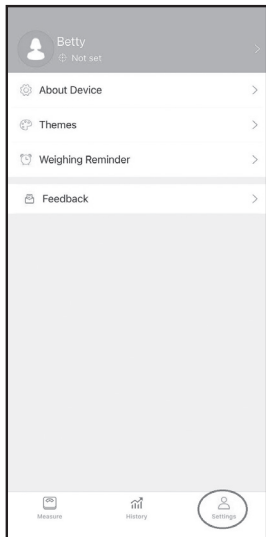Due to the continuous product improvement and updates on the Mobile App interface, the above mention data and user interfaces are only for reference. The actual APP interface is subjected to the final download version and won't notice to user in advance.

## OTHER OPERATION INSTRUCTIONS

- “C” displayed on the screen indicates that the weight is between 0.2-3.0kgs, information cannot be stored and the device will be reset to zero.
- Overload indication is “Err” on screen.
- Low Battery indication is “Lo” on screen.

## TECHNICAL DATA

Item no: GFS90S

Power: 3 x AAA batteries (Included)

Product Dimension: 25(H) x 280(W) x 280(D)mm

Maximum Capacity: 396 lb / 180 kg

\*The final product specification may be varied according to a different lot of material and actual input voltage power of the usage area.

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