简体中文	1
繁体中文	6
English	11
日本語	18
한국어	23
Polskie	29
Deutsche	35
Русский	41
Français	47
Português	53
Español	59

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10

Ocean Sonic Electric Toothbrush User Manual


Thanks for using Ocean Sonic Electric Toothbrush.

Important Safety Notice

- When using the product for the first time, due to the high-speed swinging of brush, if it feels stimulating or numb, please do not worry, it is normal and the feeling will weaken and disappear gradually in around 7-14 days.
- If not using for long time (the product battery has no power due to the standby power consumption), when charging for the toothbrush again, the charging indicator doesn't show, please wait for 6-10 minutes patiently, it can recover normal charging.
- If the main part of toothbrush becomes slightly hot when charging, it is normal. If you feel hot with your hand, please immediately stop charging and submit the toothbrush to Ocean's After-sales Services Center for testing.
- The product can only use the special brush head for Ocean Sonic electric toothbrush. Ocean will not be held liable for all the consequences caused by the use of other unofficial or non-licensed brush head products.
- The product is only used for cleaning teeth and must not be used for other purposes.
- The product is only intended as a personal cleaning product. It is not recommended to use it by multiple persons simultaneously.
- The product is only used indoors, do not place the toothbrush under sun blazing sun for quite a long time. It shall be placed in a dry and well-ventilated place with normal temperature.

Attention should be paid for those with the following conditions of oral sensitivities

- Those who have recently undergone periodontal surgery or are undergoing treatments.
- Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy bleeding from the gums.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with the above group or consult a professional dentist before use.
- This product is not suitable for people with disabilities, mental illness or lack of relevant experience and knowledge (including children (0-12 years old)). Unless used under the supervision of their guardians.

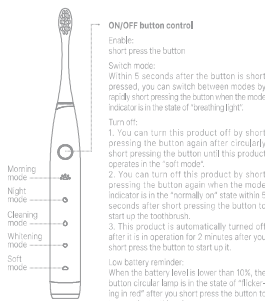11

Instruction for use

Default Brushing program

- Five toothbrushing modes are built in this product: Morning mode, night mode, cleaning mode, whitening mode and soft mode.
- Morning mode: compared with "night mode", the pressure is low. As a cleaning and care tool, it can effectively remove dental plaque caused during the sleep.
- Night mode: deeply remove food residues and dental plaque on the teeth with high pressure.
- Cleaning mode: meet requirements for routine and ordinary cleaning. Compared with "night mode", the pressure is relatively tender.
- Whitening mode: effectively clean bite marks in those who have habits of smoking or frequently drink coffee, tea, etc.
- Soft mode: protect the health of gingiva while fulfilling the fundamental function of removing dental plaque in those who have oral allergy.

Type-C direct connection charging

Open the bottom cover of the electric toothbrush, connect Type-C interface into the port at the bottom of the toothbrush, connect one end of the USB interface to the power supply, the button circular lamp flickers in white and this indicates that the toothbrush is being charged.

Charging process

Charging status: the button circular lamp is on in white and flickers.

Use for the first time


The user may charge the toothbrush depending on his or her need for battery life to ensure normal use of this product.

Application process

When the battery level is lower than 10%, the button circular lamp is in the state of "flashing in red" after you short press the button to turn on the toothbrush.

Precautions

Battery can not be replaced by users themselves. It can only be replaced by the personnel in the customer center authorized by Ocean, non-specialized personnel can not dismantle the product by themselves, or it will be regarded as that the user give up the warranty automatically. Please charge the Ocean toothbrush with a 5V/1A or 5V/2A mobile phone charger. Such mobile phone charger includes iPhone charger, tablet charger etc. If the charging current is lower than 1A, this may cause the charging time to be extended; if you charge the toothbrush with a mobile phone charger with a voltage of higher than 5V, this may cause the charger circuit to be damaged.

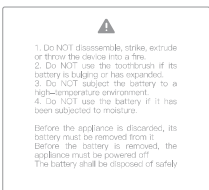12

Recommended tooth brushing method

Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.

Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.

Please gargle with clear water after brushing.

Rinse with tap water and air dry, and do not pull the bristles of the brush head

Product specifications

Product Name	Ocean Sonic Electric Toothbrush
Model	F5002
F Rating	IPX7
Power of Subject	5W
Charger	Type-C direct Recharge
Input	5V ≙ 1A
Service Time*	≙ 180 Days(*)
Battery Capacity	2500 mAh(**)
Charging time	≙ 5 hours(***)
Implementation standard	GB 4706.1-2005 GB 4706.69-2008

- * Validated by 2 brushing processes per day and 2 minutes of each process (Morning mode)
- ** The battery capacity theoretically has a capacity tolerance of ±5.0mAh.
- *** It is suggested to adopt the USB charger of 5V ≙ 1A in the phone industry.

13

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 -- Reorient or relocate the receiving antenna.
 -- Increase the separation between the equipment and receiver.
 -- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 -- Consult the dealer or an experienced radio/TV technician for help.
 The device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure condition without restriction.

IC WARNING

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:
 (1) This device may not cause interference.
 (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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CAN ICES-5(B)/NMB-5(B)

Hereby, Shenzhen Yunding Information Technology Co., Ltd. declares that the radio equipment type F5002 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://www.ocean.com/en/services.aspx>

14

Please request a maintenance report from the maintenance center and keep it together with the warranty certificate.

Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers a guarantee for consumers who purchase the Ocean Sonic Electric Toothbrush. In the case of any inconsistency between the service articles and the provisions of guarantee, the provisions of the guarantee and the contents of the guarantee certificate shall govern.

If the product you have purchased is found to have a fault or defect, please contact our service line at service@ocean.com

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.
 Manufacturer's address: #28G, Building 3, Dezhong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China
 Production address: Building A, No.6 of Xinnu Road, Xinnu Community, Pinghu, Shenzhen, China
 Website: www.ocean.com

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16

Warranty

- This product has the guarantee of toxic and hazardous substance or element in compliance with the requirements of the European Union in "RoHS" (2002/95/EC).
- This product has the guarantee of energy efficiency in compliance with the requirements of the European Union in "Energy Label" (2002/96/EC).
- * Except for any, four materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence as long as the environmental protection period.

Guarantee Certificate of Ocean Sonic Electric Toothbrush

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